

Support document

Example of a purification system using MassivEV™ EV Purification Column PS

Purification System (Sample Reaction)



Peristaltic pump

e.g. Repligen, Product Number: ACJR-U10-R

Sample

Cell culture supernatant

Fluid transfer tubing

Connection examples are shown at the bottom of the page


MassivEV™ EV Purification Column PS

FUJIFILM Wako, Product Number: 131-19491 / 137-19493

Collection tube

e.g. Corning, Product Number: 430791 (15 mL)

Connection examples of fluid transfer tubes

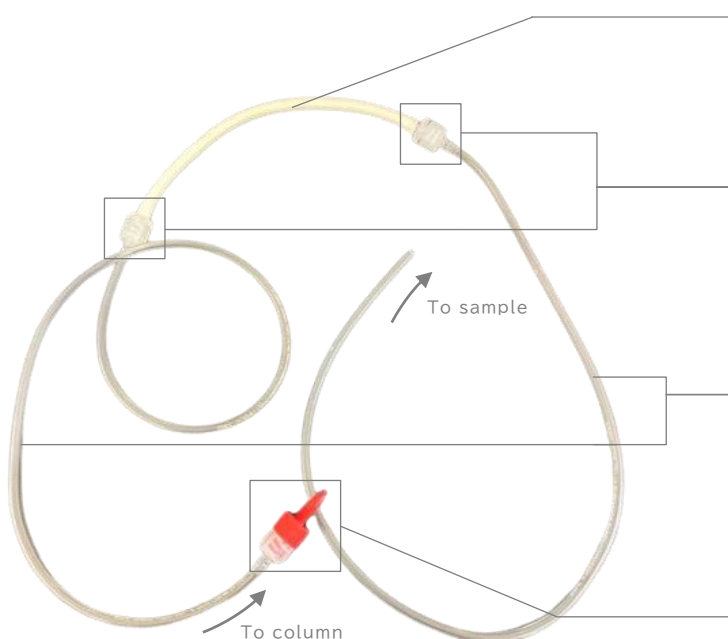

PharmaPure® tubing

e.g. Product Number: 06435-14

* Use the product specified by the manufacturer of the peristaltic pump because it affects flow velocity control.

Connection barb

e.g. Bio-Rad, Product Number: 7318222 / 7318225

TYGON® tubing

e.g. Bio-Rad, Product Number: 7318215

Connector

e.g. Cytiva, Product Number: 18111251

PharmaPure® and TYGON® are registered trademark of Saint-Gobain Performance Plastics Corporation.

This information is based on our experiments and does not guarantee performance.

Fujifilm Wako does not assume any responsibility or liability for any damage, whether direct or indirect, relating to the use of this publication.